

## PRODUCT DATA SHEET

Standing: 2024-12-02

LAUDA Variocool VC 3000

Process thermostat 230 V; 50 Hz

Part Number: L000947

### Features

- Process thermostat suitable for use with non-flammable heat transfer liquids
- Coloured TFT display for simultaneous indication of actual & set values and graphic illustration of the temperature profile
- Clear text menu navigation, six selectable languages DE, EN, FR, ES, IT, RU
- Easy input via cursor and soft keys
- Fully electronic continuous controller with PID action
- Electronic level indication and low level alarm
- Powerful pressure pump
- USB interface as standard
- Remote fault indication through floating contact
- Upgradeable with an interface module (analogue module, contact module, RS 232/485 module, Profibus, Ethernet-USB module)
- Integrated programmer with max. 150 segments, splittable in 5 programmes
- Adjustable bypass for pressure limiting
- Filler opening on top, drain tap on the backside
- SmartCool system for energy-saving digital cooling management including compressor on-off control
- Utilises traditional refrigerants (HFCs) in accordance with European legislation to control F-gases (EU) 573/2024
- Condenser cooling Air
- Operates with non flammable liquids (water, water/glycol)



Reserve technical changes



Working temperature min.  
-20 °C



Working temperature max.  
80 °C

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser

# PRODUCT DATA SHEET

Standing: 2024-12-02

LAUDA Variocool VC 3000

Process thermostat 230 V; 50 Hz

Part Number: L000947

## Technical Features (according to DIN 12876)

Working temperature range	-20 ... 80 °C
Ambient temperature range	5 ... 40 °C
Temperature stability	0.05 ± K
Heater power max.	1.5 kW
Power consumption max.	2.6 kW
Current max.	14 A
Pump Pressure max.	4,8 bar
Pump flow rate max. (pressure)	37 L/min
In / Outlet connection thread (outside)	G 3/4"
Pressure adjustment	bypass
Filling volume max.	33 L
Overall dimensions (WxDxH)	550 x 650 x 970 mm
Weight	89 kg
Refrigerant stage 1	R-449A (GWP 1397); 0.950 kg; 1.3 t CO <sub>2</sub> -eq
Power supply	230 V; 50 Hz
Power plug	Power cord with plug (SEV 1011, SEV 5934/2, T23)

Reserve technical changes

Temperature	Heat transfer liquid	Cooling Capacity 50 Hz
20 °C	Ethanol	2.8 kW
10 °C	Ethanol	2.2 kW
0 °C	Ethanol	1.48 kW
-10 °C	Ethanol	0.75 kW
-20 °C	Ethanol	0.25 kW

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser

## PRODUCT DATA SHEET

Standing: 2024-12-02

LAUDA Variocool VC 3000

Process thermostat 230 V; 50 Hz

Part Number: L000947

### Standard accessories

- 2 nipples 3/4" with screw cap G3/4 for pump connectors

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser

## PRODUCT DATA SHEET

Standing: 2024-12-02


LAUDA Variocool VC 3000

Process thermostat 230 V; 50 Hz

Part Number: L000947

### PUMP CHARACTERISTIC Heat transfer liquid: Water

Pressure bar

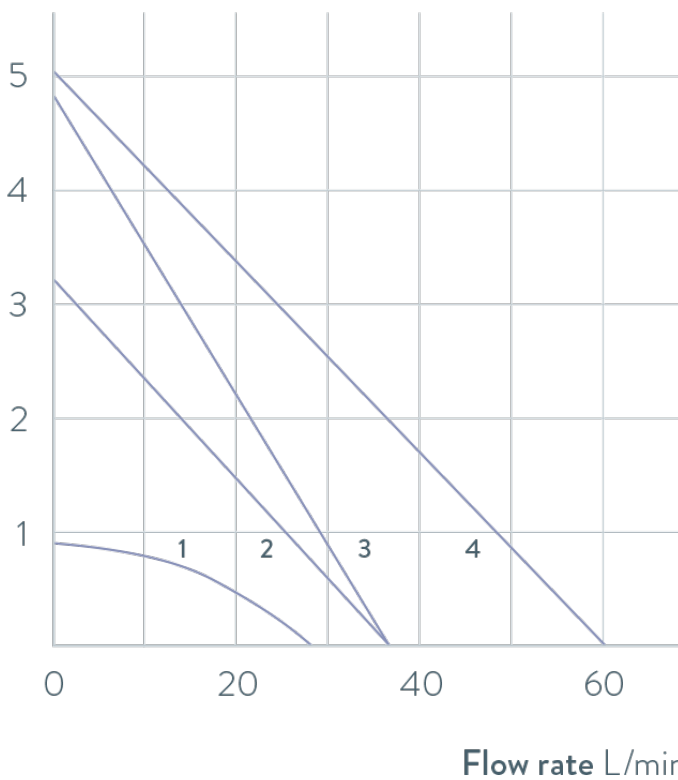

- 1 0.9 bar, 28 L/min
- 2 3.2 bar, 37 L/min
- 3 4.8 bar, 37 L/min
- 4 5.0 bar, 60 L/min

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser